




Available in Roth® Type*, MBT™*, .020" x .028" and Extra Torque.
Packaged in single patient kits, sliding tubes for A-Box or in bags of 5.

Visit www.gcorthodontics.eu for further information.



Printed on a 60% recycled paper. - 99-11E0-1100 Rev01 - 12/2015 * EXPERIENCE is a trademark of GC Orthodontics Europe GmbH. All other trademarks are property of their respective holders. All Rights Reserved.

EXPERIENCE ceramic

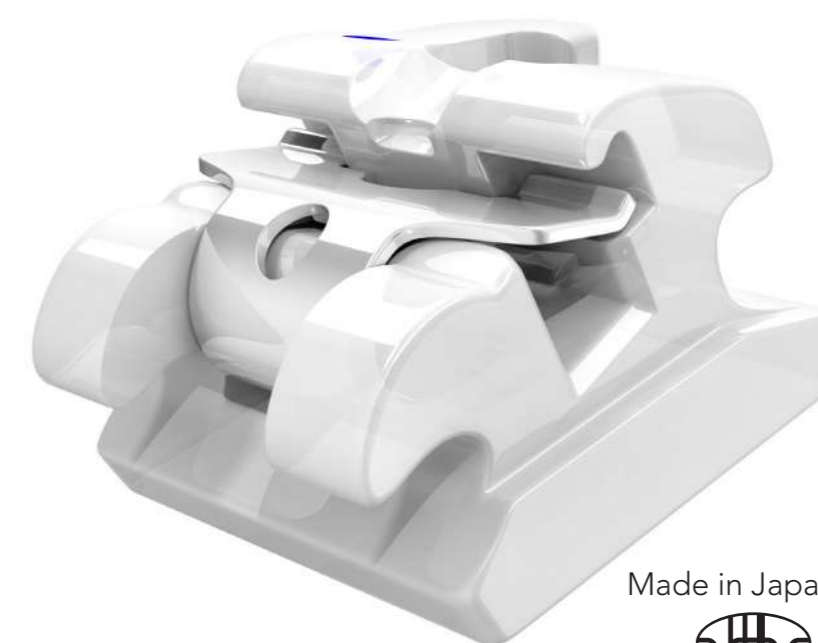

Made in Japan by



Reinvented with a **New Clip**

GC Orthodontics
GLOBAL EXPERTS IN ORTHODONTICS

EXPERIENCE ceramic

SL Product **Hi Aesthetic**


Discreet when in the mouth, this ceramic injected molded self-ligating aesthetic appliance features all the mechanical properties and strength of metal brackets. Highly translucent, its rhodium coated clip literally blends into the bracket's ceramic body.

Rhodium Coated **Rh**

Rhodium coated technology is also used in our Initialloy and BIO-ACTIVE archwires. When combined with EXPERIENCE™ ceramic brackets, the overall aesthetics of the appliance are exceptional.

Color Codes

Disto-gingival color coding allows for quick bracket identification and bonding accuracy. Color codes may differ for bicuspid for some prescriptions.

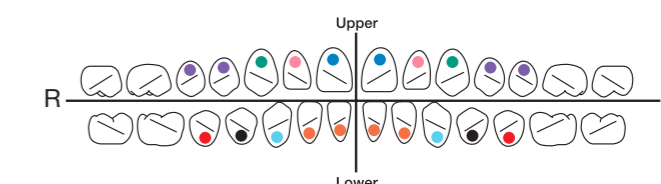

Gingival clip-opening system

Disto-gingival color ID

Interactive clip

Triple chamfered slot entry.
Milled Slot for an incredible precision

Molded mechanical base
for 100% safe retention
and easy debonding

New hook design for cuspid
and premolar brackets

Buccal clip-opening system

Spring clip with full slot
coverage and two contact
points

C-Tool: Features both a double forked
end to facilitate wire insertion and an
end with a hollow ball to safely open the
bracket from the gingival side.

P/N. 50-1130-0000

Open Tool: Is a double-ended instrument
that features a hollow ball, which aides in
opening EXPERIENCE brackets from either
the buccal (hook) or gingival (safety) side.

P/N. 50-0030-0000